




# Shakers/Mixers/Rollers

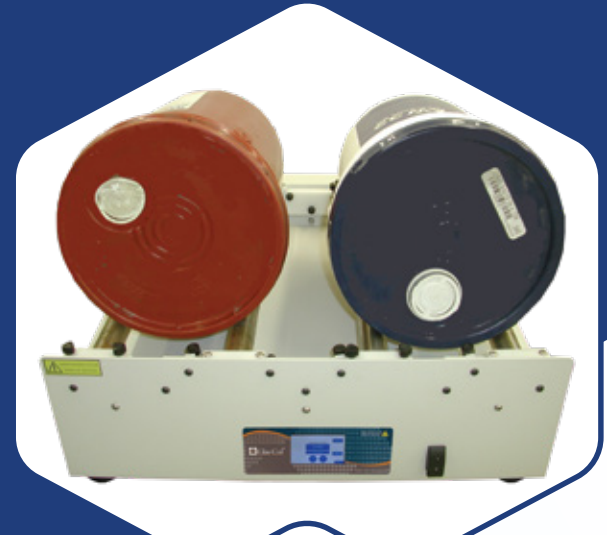


Catalog Number	Voltage	Description
099A 0616700244	120V	Digital - 27"W x 31"D x 10"H Roller
099A 0616700424	120V	Analog - 27"W x 31"D x 10"H Roller
099A 0616700357	240V	Digital - 27"W x 31"D x 10"H Roller

## Roller:

The Glas-Col roller is designed for laboratory or industrial applications. The heavy-duty construction accommodates up to two 5-gallon industrial containers, six 1-gallon cans, or 8-quart can/jars. In addition it can roll any size bottle 4-1/4" diameter and larger. The Digital version comes with a touch screen interface and the Analog version a speed adjust knob.

There is an optional cabinet that can be designed to contain the roller. A cutout in the front to gives the operator easy access the control panel on the roller. A door on the top of the cabinet opens for the loading and unloading of the containers.

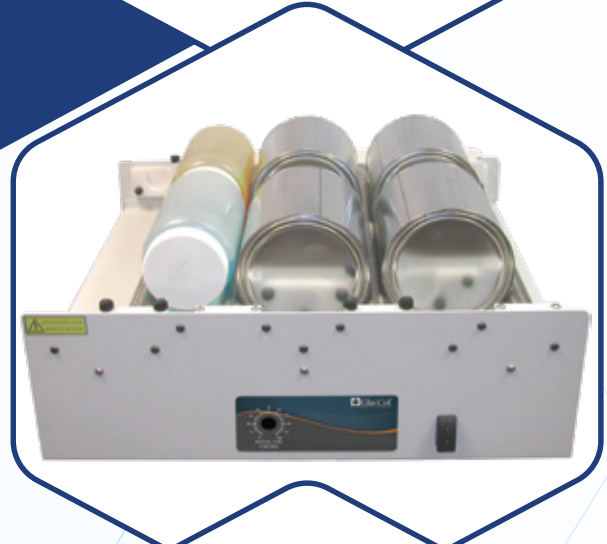

## Specifications:

- Digital system: Touchscreen interface for control of speed and time functions.
- Several alarm conditions available.
- Create up to 10 recipes.
- Analog system: Speed adjust knob.

## Typical speed of containers:

5 gallon:	20 rpm max
2 gallon:	55 rpm max
1 gallon:	38 rpm max
Quart:	58 rpm max

